

# Gum Hero™

Perio Home Care System



Gum Hero™ is a personalized perio home care system designed to fight harmful bacteria in the mouth that causes gum disease.

The Gum Hero™ Perio Home Care System consists of:

- ▶ Personalized trays designed for each patient, ensuring a precise fit around the teeth and gums.
- ▶ An oral cleansing gel which is applied inside the tray to fight bacteria below the gumline, where brushing and flossing can't reach.




## Periodontal Pocket Depth

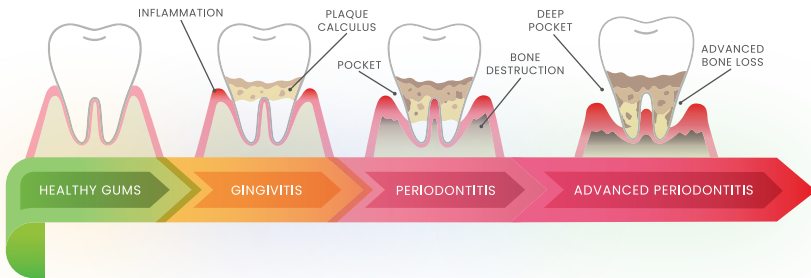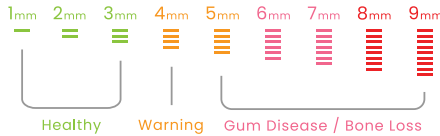

Your dentist or hygienist may use a tool called a periodontal probe to measure the spaces between your teeth and gums. These spaces are called **periodontal pockets**.

Deeper pockets may be a sign of uncontrolled gum disease. Fortunately, the Gum Hero™ perio home care system can help deliver oral cleansing gel into these areas to fight harmful bacteria.

If left uncontrolled, gum disease can lead to gum redness, inflammation, bleeding, and tooth loss.

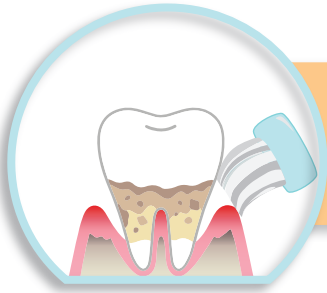

Gum Hero™ is a perio home care system designed to fight harmful bacteria deep below the gumline where brushing and flossing can't reach.

## How to use Gum Hero™ trays:

- ▶ Apply a thin line of Gum Hero™ gel to the inside of the trays
- ▶ Firmly place the trays over your teeth and gently bite down to seat them. Spit out any excess gel
- ▶ Wear the trays for a minimum of 15 minutes, or as directed by your dentist or healthcare provider
- ▶ Remove the trays and rinse your mouth to remove any excess gel
- ▶ Clean the trays with cold water and place them in the provided tray caddy to dry



Daily usage of the Gum Hero™ perio home care system will help fight harmful bacteria in the mouth and control the symptoms of gum disease such as redness, inflammation, and bleeding.

Healthy Gum

Diseased Gum

